


N-Glycan / Labeled N-Glycan based on Chemical Synthesis


N-Glycan

Our chemical synthesis of oligosaccharide using various sugar building blocks provides various structure-defined *N*-glycans such as human type glycans (G0, G2 and SG), uniform isomer (each G1 and G1') and heterogenic antigen (including α Gal or NeuGc).

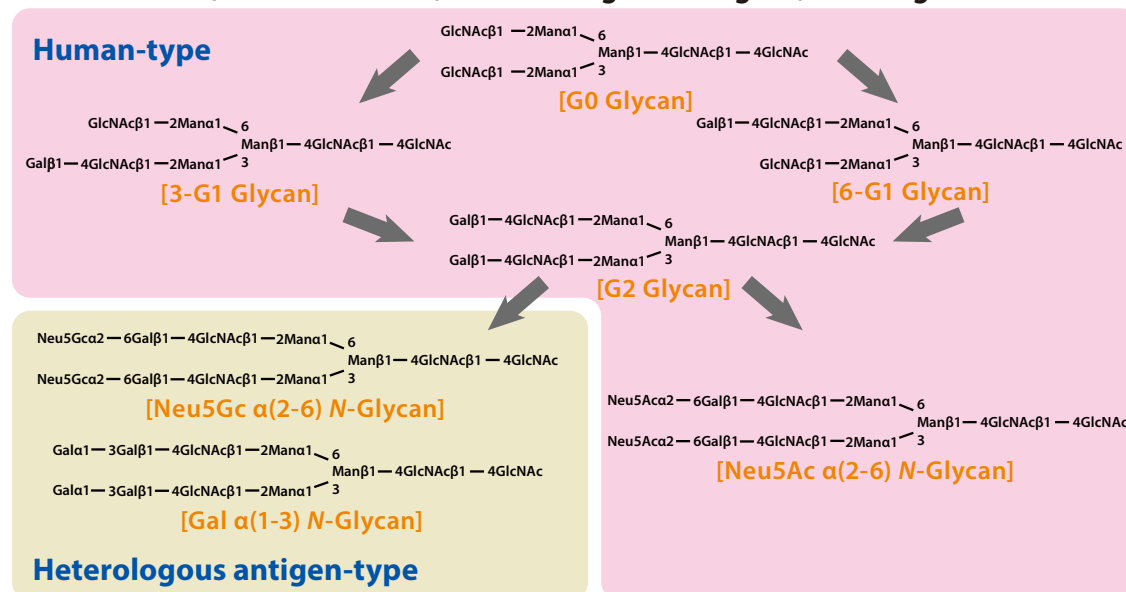


Labeled N-Glycan

We produce high grade 2-AB labeled *N*-glycans. Please take advantage of our fluorescent-labeled product portfolio as a test sample for MS, CE and HPLC analyses.

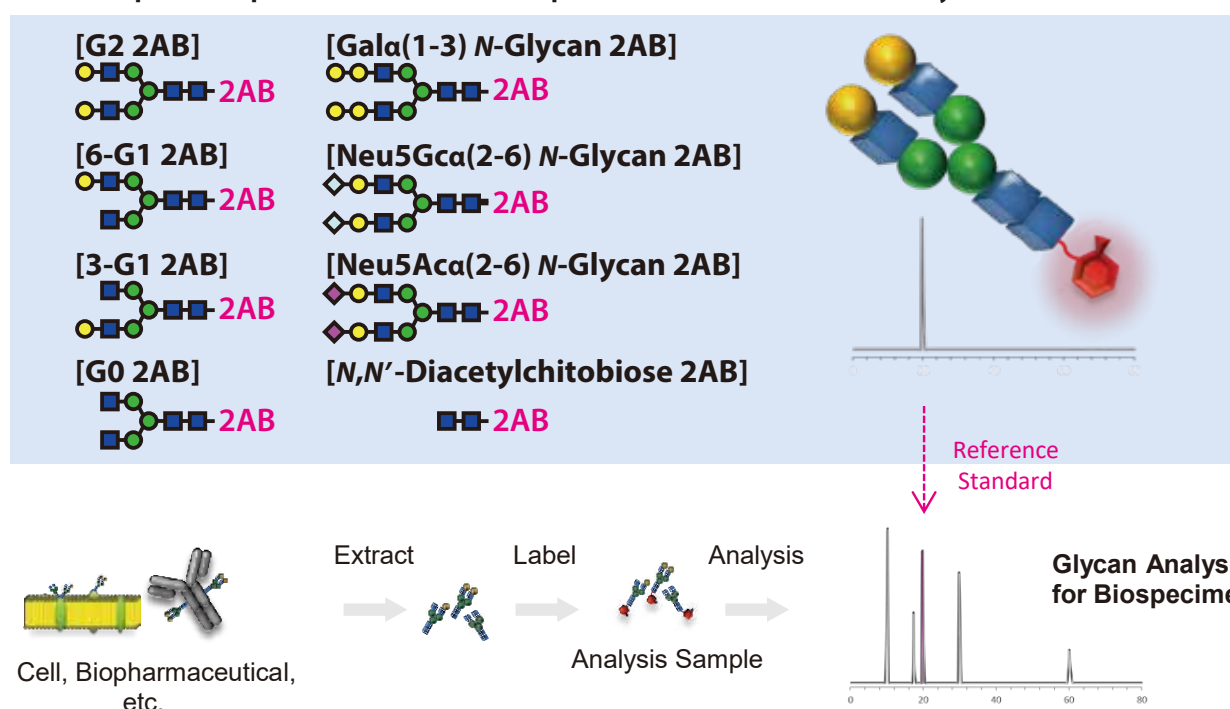

Application of 2AB-labeled Glycans

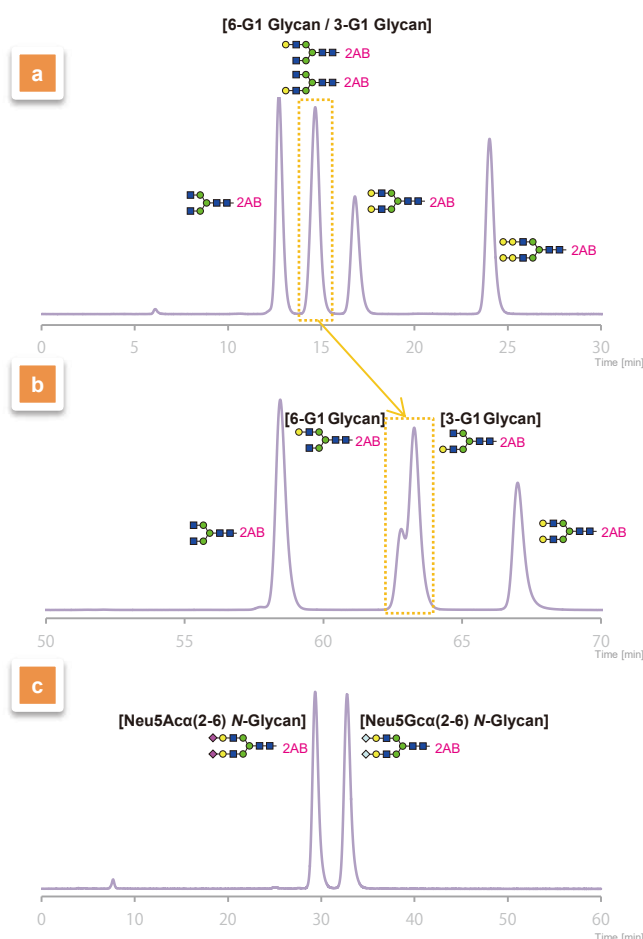

[Column] Asahipak NH2P-50 4E (Φ4.6×250 mm)

[Detection] FL (Ex: 330 nm, Em: 420 nm)

[Inject] 5 μ L (10 μ g/mL except for 5 μ g/mL of 6-G1 2AB)

[Flow rate] 0.5 mL/min

[Temp.] 40°C

[Solvent] A : Acetonitrile

B : 50 mM Ammonium formate (pH4.4)

[Condition]

- (A) 0-10 min B: 40%
10-30 min B: 40% (gradient)→50%
- (B) 0-20 min B: 20%
20-60 min B: 20% (gradient)→40%
60-80 min B: 20%
- (C) 0-60 min B: 100%

(A) The labeled neutral *N*-glycan mixture was individually separated based on the distinctive features of glycan structure. (However, each G1 isomer was not separated.)

(B) Each G1 isomer; 6-G1 and 3-G1, is modestly detached in this condition.


(C) The acidic *N*-glycan; Heterologous antigen-type Neu5Gc α (2-6) *N*-Glycan was completely separated from Neu5Aca(2-6) *N*-Glycan due to the minimal difference of sialic acid.

N-Glycan / Labeled N-Glycan based on Chemical Synthesis

G0 Glycan	[G0484]	G0 2AB	[G0490]
3-G1 Glycan	[G0485]	3-G1 2AB	[G0491]
6-G1 Glycan	[G0486]	6-G1 2AB	[G0492]
G2 Glycan	[G0487]	G2 2AB	[G0493]
Neu5Aca(2-6) <i>N</i> -Glycan	[N1065]	Neu5Aca(2-6) <i>N</i> -Glycan 2AB	[N1073]
Neu5Gca(2-6) <i>N</i> -Glycan	[N1064]	Neu5Gca(2-6) <i>N</i> -Glycan 2AB	[N1075]
Gal α (1-3) <i>N</i> -Glycan	[G0488]	Gal α (1-3) <i>N</i> -Glycan 2AB	[G0494]
		<i>N,N'</i> -Diacetylchitobiose 2AB	[D5327]

A part of the products was commercialized by support from the projects of New Energy and Industrial Technology Development Organization (NEDO) and the Japan Agency for Medical Research and Development (AMED) implementation of manufacturing technology association of biologics.


N-Glycan / Labeled N-Glycan

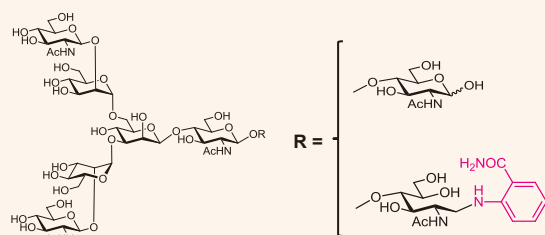

G0 Glycan [G0484]




G0 2AB [G0490]

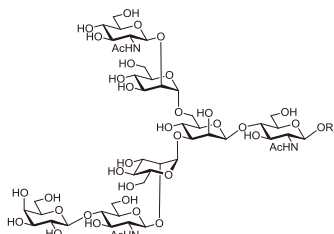

3-G1 Glycan [G0485]




3-G1 2AB [G0491]

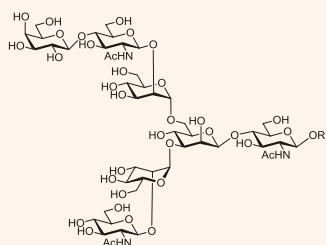

6-G1 Glycan [G0486]




6-G1 2AB [G0492]

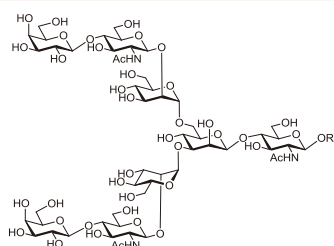

G2 Glycan [G0487]




G2 2AB [G0493]

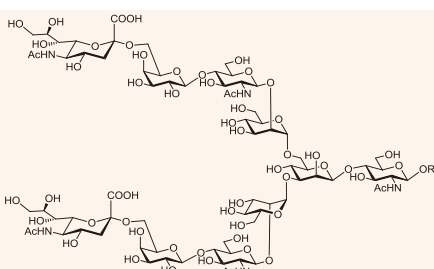

Neu5Aca(2-6) N-Glycan [N1065]




Neu5Aca(2-6) N-Glycan 2AB [N1073]

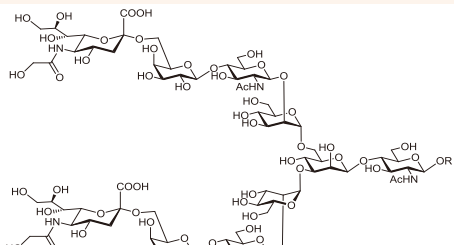

Neu5Gca(2-6) N-Glycan [N1064]



Neu5Gca(2-6) N-Glycan 2AB [N1075]

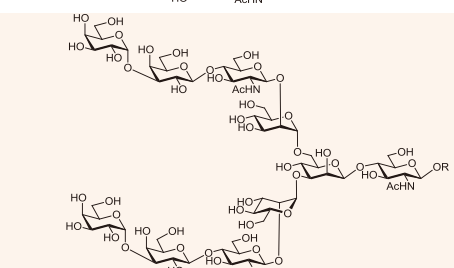

Galα(1-3) N-Glycan [G0488]





Galα(1-3) N-Glycan 2AB [G0494]





NMR Data of Labeled *N*-Glycan

¹H NMR JEOL Ltd.
[Solvent] D₂O
[Internal Standard] *tert*-BuOH



G0 2AB

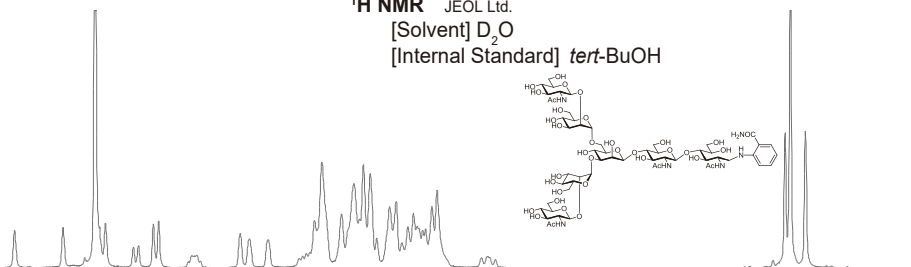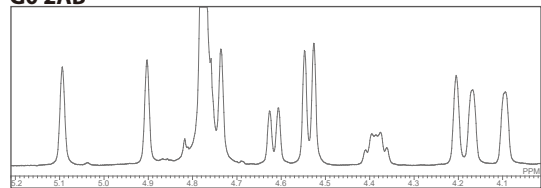



3-G1 2AB

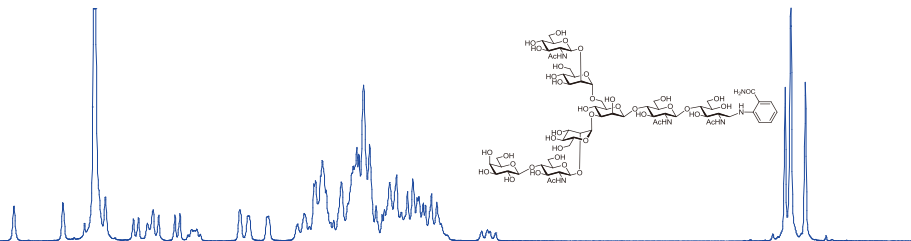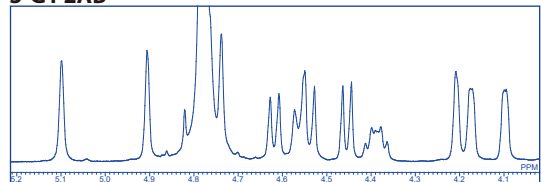



6-G1 2AB

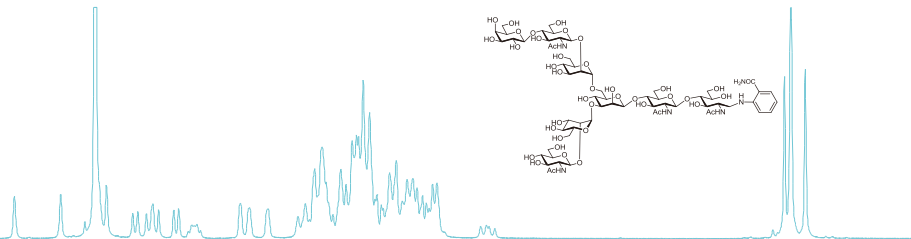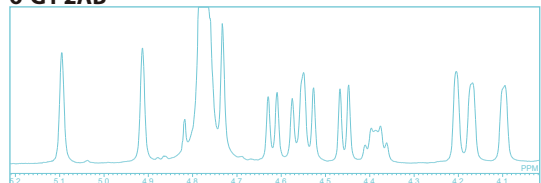



G2 2AB

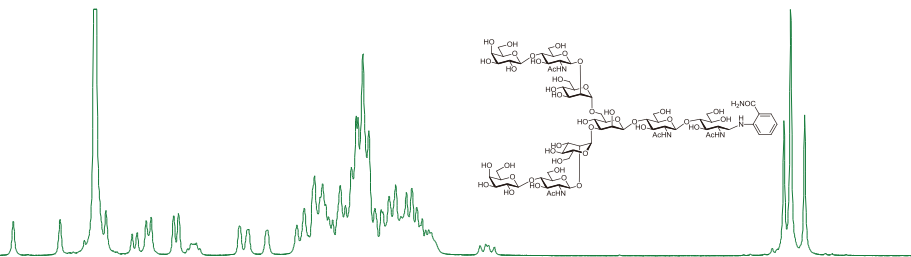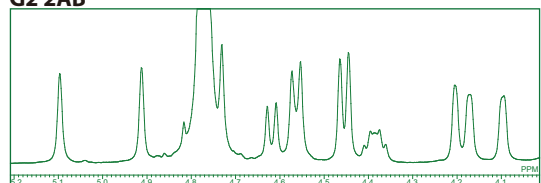



Neu5Acα(2-6) *N*-Glycan 2AB

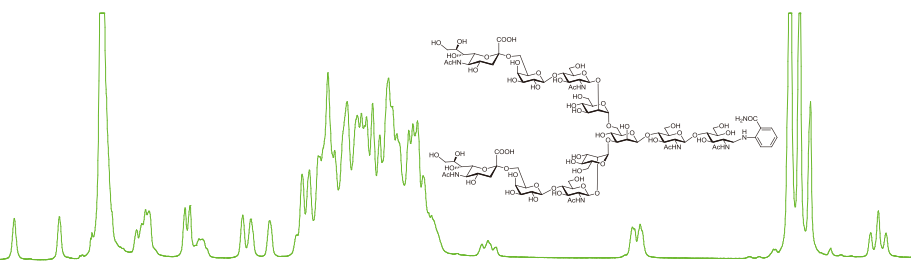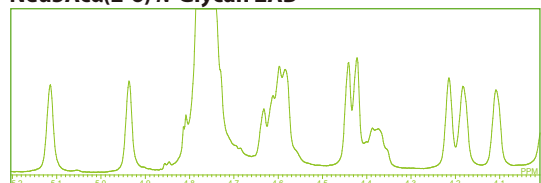

Neu5Gcα(2-6) *N*-Glycan 2AB

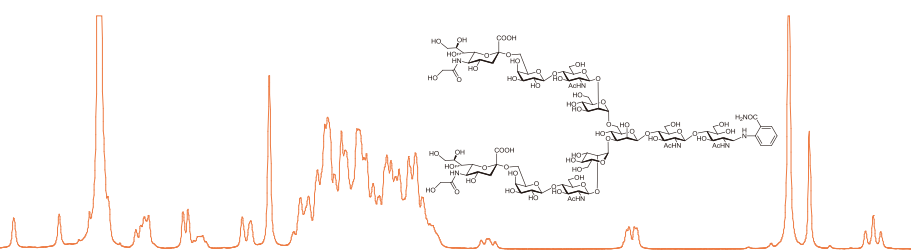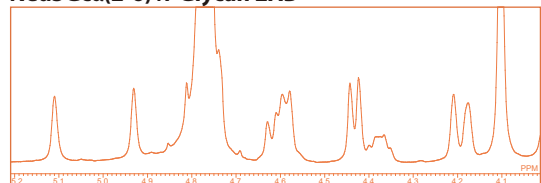

Galα(1-3) *N*-Glycan 2AB

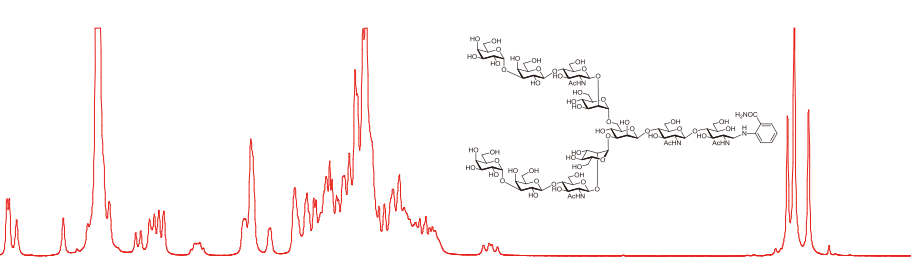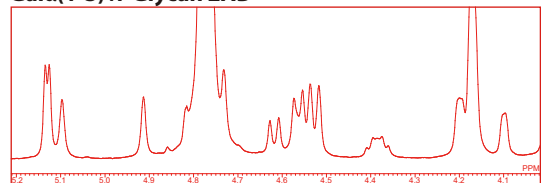

7.0

6.0

5.0

4.0

3.0

2.0








PPM

MS / HPLC Data of Labeled *N*-Glycan

MALDI-TOF Bruker Daltonics

[Mode] Positive mode

(Neu5Acq(2-6) *N*-Glycan 2AB, Neu5Gcq(2-6) *N*-Glycan 2AB were analyzed by Negative mode)

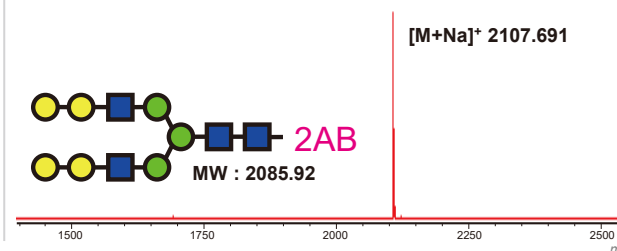
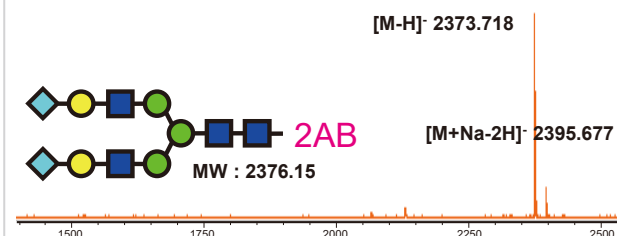
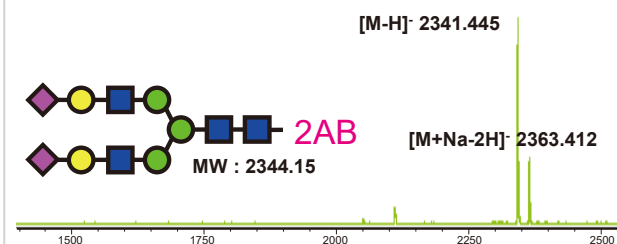
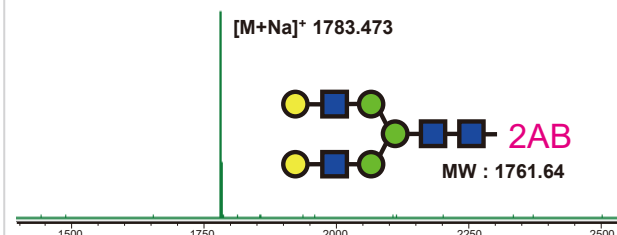
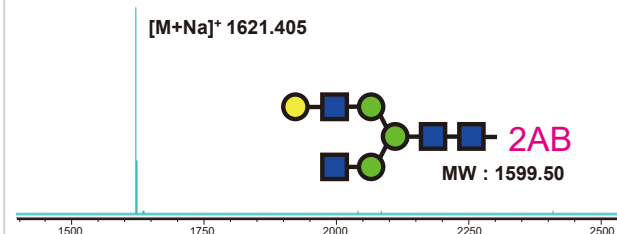
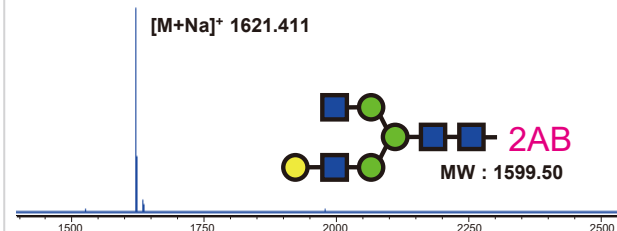
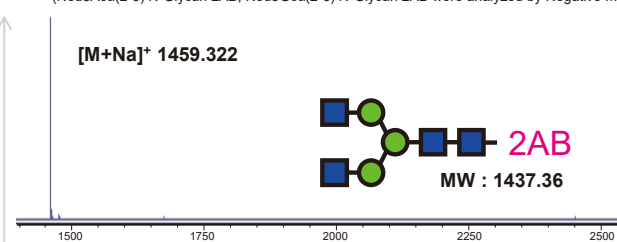


HPLC Hitachi Hightech

[Detection] FL (Ex: 330 nm, Em: 420 nm)

[Column] Asahipak NH2P-50 4E (4.6×250 mm)


HPLC purity >95%

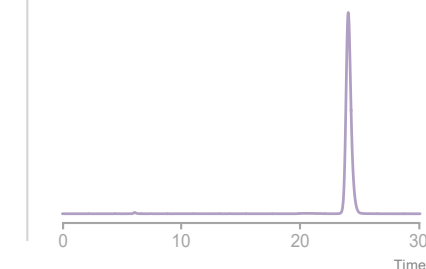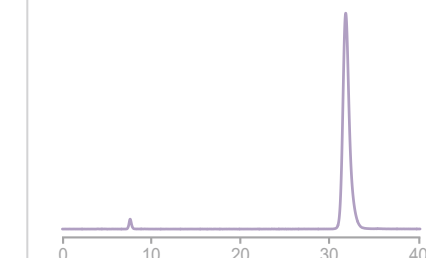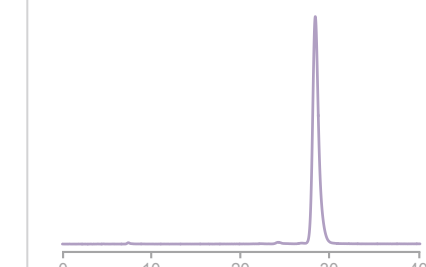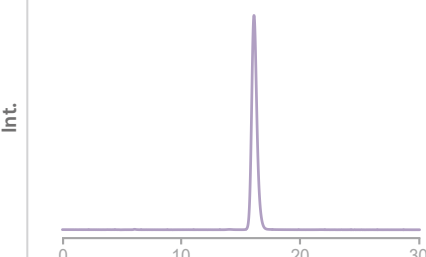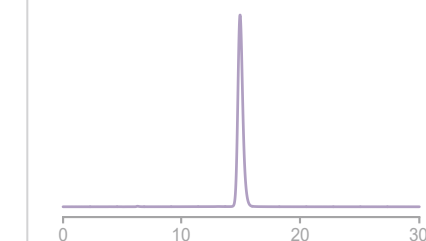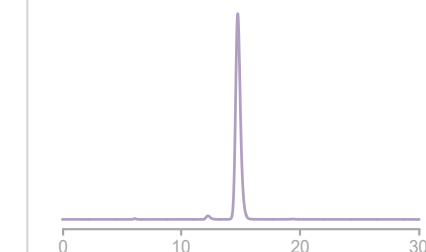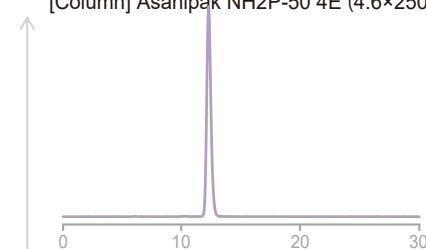

Human-type

Heterologous antigen-type

Glycoscience

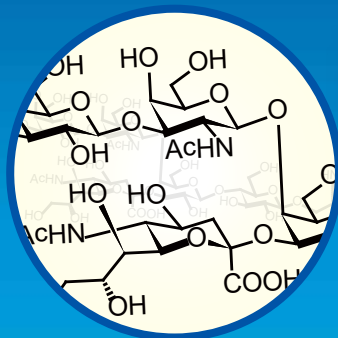


Sugar chains

- Functional oligosaccharides
- Various sugar blocks at 10-100kg
- Synthetic technology with high-quality
- Application of sugar-conjugates



Enzymes

- Endo type enzymes for glycoconjugate synthesis
- Glycohydrolase and their substrates

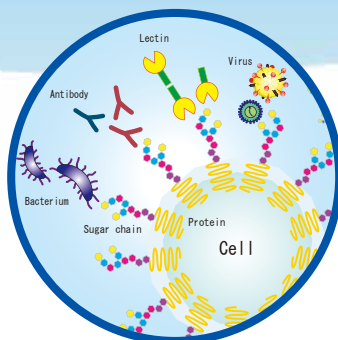

Antibodies & Lectins

- Antibodies for glycolipids and glycosaminoglycans
- Fucose specific lectin
- Modification of antibodies and lectins

Ordering and Customer Service

TCI AMERICA

Tel : 800-423-8616 / 503-283-1681
Fax : 888-520-1075 / 503-283-1987
E-mail : Sales-US@TCIchemicals.com

TCI EUROPE N.V.

Tel : +32 (0)3 735 07 00
Fax : +32 (0)3 735 07 01
E-mail : Sales-EU@TCIchemicals.com

TCI Deutschland GmbH

Tel : +49 (0)6196 64053-00
Fax : +49 (0)6196 64053-01
E-mail : Sales-DE@TCIchemicals.com

Tokyo Chemical Industry UK Ltd.

Tel : +44 (0)1865 78 45 60
E-mail : Sales-UK@TCIchemicals.com

梯希爱(上海)化成工业发展有限公司

Tel : 800-988-0390 / 021-67121386
Fax : 021-6712-1385
E-mail : Sales-CN@TCIchemicals.com

Tokyo Chemical Industry (India) Pvt. Ltd.

Tel : 1800 425 7889 / 044-2262 0909
E-mail : Sales-IN@TCIchemicals.com

TOKYO CHEMICAL INDUSTRY CO., LTD.

Tel : +81 (0)3-5640-8878
E-mail : globalbusiness@TCIchemicals.com

• Chemicals itemized in this brochure are for research and testing use only. Please avoid use other than by chemically knowledgeable professionals. • Information such as listed products and its specifications and so on are subject to change without prior notice. • The contents may not be reproduced or duplicated in whole or in part without permission of Tokyo Chemical Industry Co., Ltd.