


# High-purity Silane Coupling Agent: Octadecyltrichlorosilane (>99.0%)

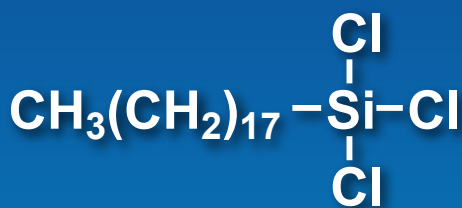



[T3815]

### Advantages

- Extremely low levels of any branched and short alkyl chain impurities
- Packaged in a sealed ampule under argon to maintain high purity
- High-quality surface treatment enhances device performances

## Comparison of OFET performance by purity of octadecyltrichlorosilane (ODTS)

This product [T3815] and the conventional product [O0079] were evaluated as surface treatment agents in organic field effect transistor (OFET) devices using an in-house system.

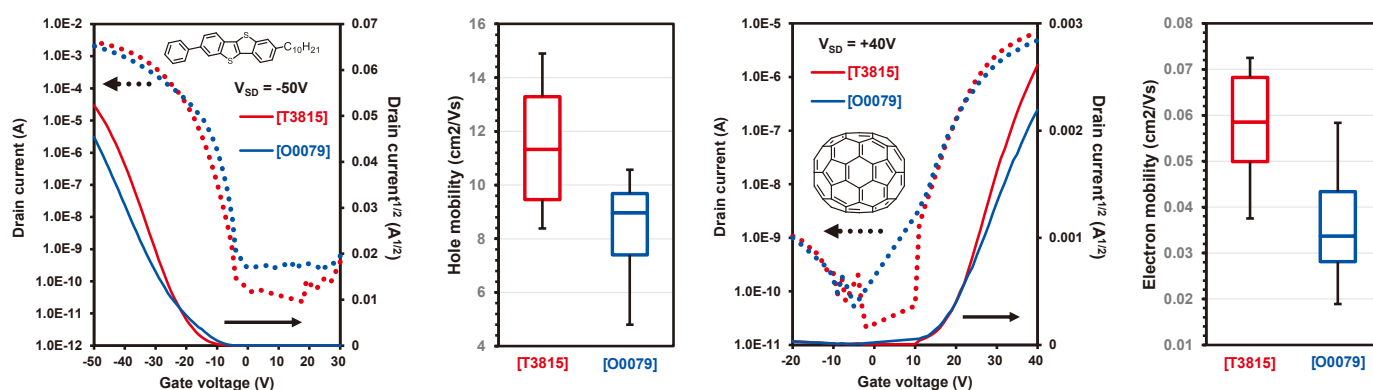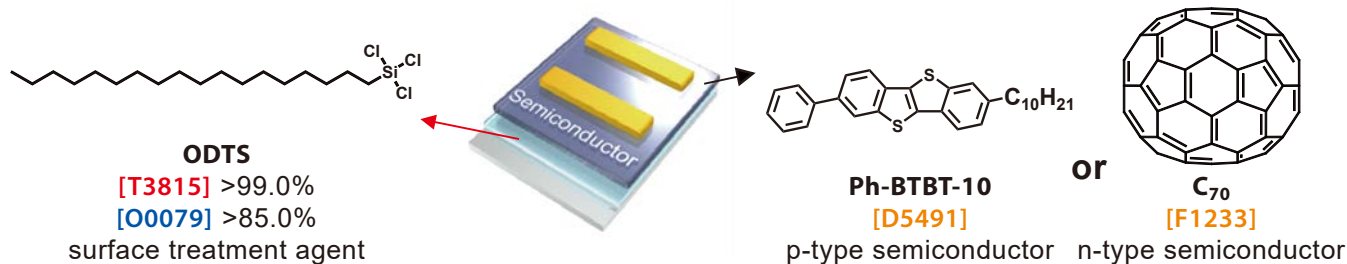

	ODTS		$\mu$ ave.	$\mu$ max	$V_{Th}$	on/off
			[cm <sup>2</sup> /Vs]	[cm <sup>2</sup> /Vs]	[V]	
Ph-BTBT-10 [D5491]	[T3815]	>99.0%	11.3 (2.1)	14.9	-17.6 (3.7)	~10 <sup>8</sup>
	[O0079]	>85.0%	8.5 (1.6)	10.6	-19.6 (6.9)	~10 <sup>8</sup>
C <sub>70</sub> [F1233]	[T3815]	>99.0%	$5.81 \times 10^{-2}$ ( $1.12 \times 10^{-2}$ )	$7.25 \times 10^{-2}$	18.0 (1.5)	~10 <sup>6</sup>
	[O0079]	>85.0%	$3.53 \times 10^{-2}$ ( $1.00 \times 10^{-2}$ )	$5.83 \times 10^{-2}$	16.5 (2.2)	~10 <sup>6</sup>

The values in parentheses are standard deviations.


**The OFET performances were enhanced by using high-purity ODTS!**

# High-purity Silane Coupling Agent: Octadecyltrichlorosilane (>99.0%)

## Example: surface treatment procedure of Si/SiO<sub>2</sub> substrate

1. Piranha etching (H<sub>2</sub>SO<sub>4</sub>:H<sub>2</sub>O<sub>2</sub>=4:1, 80°C, 2h)
2. Ultrasonic cleaning (DW, Acetone, IPA, 10 minutes each)
3. Boiling cleaning (IPA, 3 minutes)
4. UV/O<sub>3</sub> cleaning (60 minutes)
5. ODTs solution immersion (0.01 M toluene, 16h, N<sub>2</sub>)
6. Ultrasonic cleaning (Toluene, Acetone, IPA, 10 minutes each)

### Contact angles of water

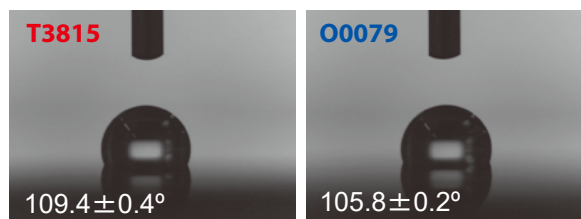

**Octadecyltrichlorosilane (= ODTs) (>99%)**

1g [T3815]

## Related products

### *n*-Alkyltrichlorosilanes

<b>Methyltrichlorosilane (C1)</b>	>98.0% (GC)(T)	25g / 100g / 500g [M0450]
<b>Ethyltrichlorosilane (C2)</b>	>98.0% (GC)(T)	25g [E0188]
<b>Propyltrichlorosilane (C3)</b>	>99.0% (GC)	25g [P0371]
<b>Butyltrichlorosilane (C4)</b>	>98.0% (GC)	25g [B0393]
<b>Pentyltrichlorosilane (C5)</b>	>97.0% (GC)	5g / 25g [A0870]
<b>Hexyltrichlorosilane (C6)</b>	>97.0% (GC)	25mL [H0547]
<b><i>n</i>-Octyltrichlorosilane (C8)</b>	>98.0% (GC)(T)	25g / 250g [O0168]
<b>Decyltrichlorosilane (C10)</b>	>98.0% (GC)	25mL [D1486]
<b>Dodecyltrichlorosilane (C12)</b>	>97.0% (T)	25mL [D1509]
<b>Hexadecyltrichlorosilane (C16)</b>	>98.0% (GC)	25g [T3532]
<b>Octadecyltrichlorosilane (C18)</b>	>85.0% (GC)	25g / 100g [O0079]

### Other-type surface treatment agent

<b>1,1,1,3,3,3-Hexamethyldisilazane (HMDS)</b>	>96.0% (GC)	25mL / 100mL / 500mL [H0089]
--	-------------	------------------------------

### Organic semiconductors

<b>Ph-BTBT-10</b>		100mg / 250mg [D5491]
<b>C<sub>70</sub></b>		100mg [F1233]

For further information please refer to our website at [www.TCIchemicals.com](http://www.TCIchemicals.com). ▶▶▶

silane coupling agent



## Ordering and Customer Service

### TCI AMERICA

Tel : 800-423-8616 / 503-283-1681  
Fax : 888-520-1075 / 503-283-1987  
E-mail : Sales-US@TCIchemicals.com

### TCI EUROPE N.V.

Tel : +32 (0)3 735 07 00  
Fax : +32 (0)3 735 07 01  
E-mail : Sales-EU@TCIchemicals.com

### TCI Deutschland GmbH

Tel : +49 (0)6196 64053-00  
Fax : +49 (0)6196 64053-01  
E-mail : Sales-DE@TCIchemicals.com

### Tokyo Chemical Industry UK Ltd.

Tel : +44 (0)1865 784560  
Fax : +44 (0)1865 784561  
E-mail : Sales-UK@TCIchemicals.com

### TCI Chemicals (India) Pvt. Ltd.

Tel : 1800 425 7889 / 044-2262 0909  
Fax : 044-2262 8902  
E-mail : Sales-IN@TCIchemicals.com

### 梯希爱(上海)化成工业发展有限公司

Tel : 800-988-0390 / 021-67121386  
Fax : 021-6712-1385  
E-mail : Sales-CN@TCIchemicals.com

### TOKYO CHEMICAL INDUSTRY CO., LTD.

Tel : +81 (0)3-5640-8878  
Fax : +81 (0)3-5640-8902  
E-mail : globalbusiness@TCIchemicals.com

Availability, price or specification of the listed products are subject to change without prior notice. Reproduction forbidden without the prior written consent of Tokyo Chemical Industry Co., Ltd.